


- Reconstructive option for large cranial defects
- Provides an alternative to customized implants, grafts, and other implant materials
- Approximates shape of the cranial hemisphere while allowing surgeon to bend and trim implant to meet the needs of the individual patient
- For easier bending and shaping, submerge the implant in a hot, sterile saline bath (82-100°C, 180-212°F)
- Screws or sutures can be placed in any desired location
- Fully porous material for tissue ingrowth
- Easy to trim with scissors or scalpel
- Individually packaged, MR Safe, Sterile, 5 Year Shelf Life



CAT#	Description	A	B	C
4116	Left	124mm	170mm	4.5mm
4117	Right	124mm	170mm	4.5mm
4118	Left	124mm	170mm	6mm
4119	Right	124mm	170mm	6mm

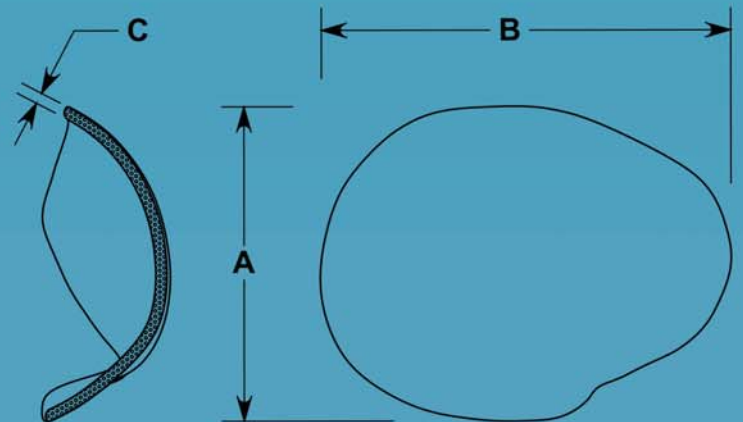