


- Reconstructive option for large cranial defects
- Provides an alternative to customized implants, grafts, and other implant materials
- Approximates contour of the superior 1/3 of the cranium
- For easier bending and shaping, submerge the implant in a hot, sterile saline bath (82-100°C, 180-212°F)
- Screws or sutures can be placed in any desired location
- Fully porous material for tissue ingrowth
- Easy to trim with scissors or scalpel
- Individually packaged, MR Safe, Sterile, 5 Year Shelf Life



CAT#	A	B	C
4120	150mm	180mm	4mm
4121	150mm	180mm	6mm

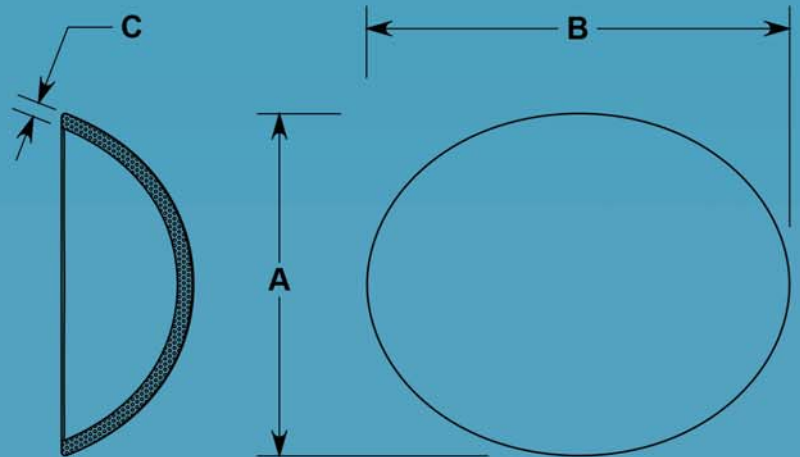