



Surgical Implants

- For repair of large, complex orbital fractures
- Design allows it to be cantilevered from the orbital rim to serve as a stable platform for orbital soft tissue when fractures are too large to support the implant in the posterior orbit
- Addition of single rigid fixation plate provides increased structural support
- Internal channel accepts microplates from conventional plating systems, facilitating fixation to bone
- Available with optional heat bonded non-porous membrane layer
- Implant bends and holds desired shape without heating
- Easy to trim with scissors or scalpel
- Individually packaged, MR Safe, Sterile, 5 Year Shelf Life

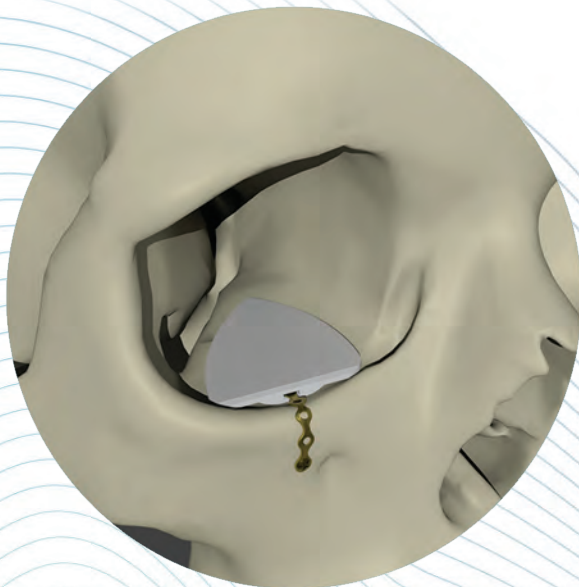

Microplate Shown. Not Included.

Shown with optional Non-Porous membrane

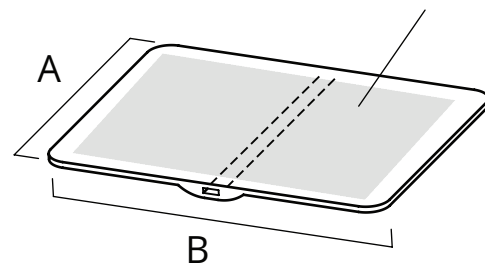

CAT#	Description	A	C
4297	Regular	38mm	50mm
4296	Membrane	38mm	50mm